



Self-retaining
Baska Bite Block
with advanced oxygenating system
for use in



BRONCHOSCOPY

Please read the instructions before use.

ADVANTAGES

- Protects bronchoscope from being bitten by the patient.
- Advanced oxygenating system.
- Prevents the Bite Block from dislodging and falling out during the procedure.
- Ease of insertion and guides the bronchoscope directly to the larynx.
- By passing the bronchoscope through the device, contamination of the tracheobronchial tree is avoided.


INSTRUCTIONS FOR USE

1. Administer 5mls of 10% Topical xylocaine syrup orally to the fasted patient to anaesthetise the oral cavity & pharynx 5 minutes before starting the procedure.
2. Place awake patient in a propped up supine position with a pillow under the head.
3. Attach 10L of O₂ to the Oxygen inlet of the stem.
4. Apply a water-soluble lubricant to the end of the stem of the plunger.
5. Insert the Bite Block between the incisors with the flange in front of the patient's teeth.
6. Preoxygenate the patient by inserting the stem halfway into the Bite Block with the O₂ connector aligned with the patient's forehead. Inserting the stem extends the Bite Block retainers so that the Bite Block remains in place during the procedure.
7. Administer IV sedation only if required just enough avoiding respiratory depression.
8. Insert the stem further far enough with O₂ port facing the forehead while avoiding contact with the posterior pharyngeal wall.
9. Then rotate the stem axially 90 degree anti-clockwise until the O₂ connector is inline with the patient's right ear. This step locks the stem in position engaging onto the ratchet.
10. Lubricate and insert the lubricated instrument through the funnel and perform procedure by advancing the flexible Bronchoscope applying local anaesthetic through the bronchoscope.
11. ETCO₂ can be monitored if required using nasal cannulas.
12. On completion of the procedure, after the instrument is retrieved, rotate O₂ connector clockwise until it is aligned with the patient's forehead. This unlocks the stem from the ratchet so that it can be withdrawn from the Bite Block.
13. Remove the stem completely to retract the Bite Block retainers which enables removal of the Bite Block with ease and oxygenate the patient with Hudson Mask until full recovery.

CAUTION

1. Single patient use only.
2. Do not use if the package is open or damaged.
3. Store in a cool dry place.
4. Avoid extreme temperatures or humidity.

Manufacturer, Distributor and Sponsor:
Baska Versatile Laryngeal Mask Pty Ltd
61 Woodside Avenue
Strathfield NSW 2135 Australia
www.baskamask.com
ABN: 42086476583

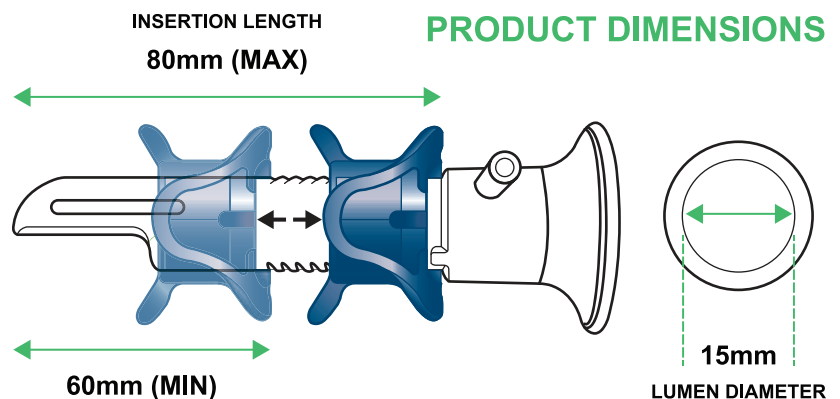

Product Order Code: B2B-B



MD



**Patents
Pending**

**ARTG ID:
468288**