



Self-retaining  
**Baska Bite Block**  
*with advanced oxygenating system*  
for use in



## GASTROSCOPY

Please read the instructions before use.

### INSTRUCTIONS FOR USE


1. Place the awake patient in the left lateral position with the head in the extended “sniffing” position.
2. Attach 10L/min O<sub>2</sub> to the oxygen inlet.
3. Insert the bite block between the incisors with the flange in front of the patient’s teeth.
4. Apply water soluble lubricant to the end of the stem.
5. Pre-oxygenate by inserting the stem halfway into the bite block with the O<sub>2</sub> connector aligned with the patient’s forehead. Inserting the stem extends the bite block retainers so that the bite block remains in place during the procedure.
6. Administer intravenous sedation and titrate the dose to maintain spontaneous respiration.
7. After the patient is adequately sedated, insert the stem further while avoiding contact with the posterior pharyngeal wall.
8. Rotate the stem 90 degrees anticlockwise until the O<sub>2</sub> connector is in line with the patient’s right ear. This locks the stem in position engaging onto the ratchet.
9. Lubricate and insert the lubricated Gastroscope through the funnel and perform the procedure.
10. ETCO<sub>2</sub> can be monitored if required using nasal cannulas.
11. On completion of endoscopy, after the scope is retrieved, rotate the O<sub>2</sub> connector clockwise until it is aligned with the patient’s forehead. This unlocks the stem from the ratchet so that it can be withdrawn from the bite block.
12. Remove the stem completely to retract the bite block retainers which enables removal of the bite block with ease and oxygenate the patient with Hudson mask until full recovery.

### CAUTION

1. Single patient use only.
2. Do not use if the package is open or damaged.
3. Store in a cool dry place.
4. Avoid extreme temperatures or humidity.

**Manufacturer, Distributor and Sponsor:**  
**Baska Versatile Laryngeal Mask Pty Ltd**  
61 Woodside Avenue  
Strathfield NSW 2135 Australia  
[www.baskamask.com](http://www.baskamask.com)  
ABN: 42086476583

### PRODUCT DIMENSIONS

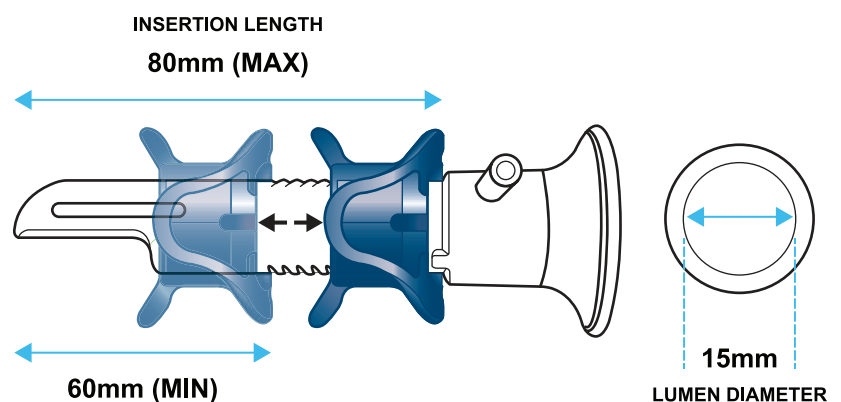

Product Order Code: B2B



**MD**



**Patents  
Pending**

**ARTG ID:  
468288**